

CERAMENT[®] G with Gentamicin
CERAMENT[®] V with Vancomycin

VALUE GUIDE

Proven bone remodeling with
local antibiotic elution

Get it right the first time[®]




THE BURDEN OF BONE INFECTION

Biofilm formation

Bone infection (osteomyelitis) is an infection of the bone that can arise from fractures, placement of implants, or diabetic foot ulcers. It can be debilitating and can lead to limb amputations.

Bacteria naturally accumulate and form sedentary colonies, commonly referred to as 'biofilm'. Once biofilm forms, susceptibility to antibiotic therapy is significantly reduced and surgical debridement is required.

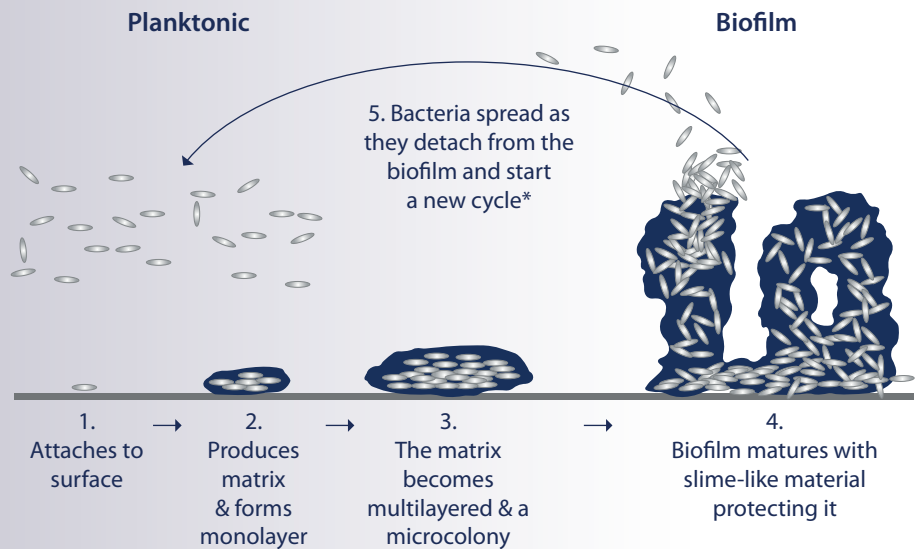

*adapted from Vasudevan et al. 2014¹

Unmet needs

Despite systemic antibiotics and medical advancements, bone infections are common and impose a significant burden:

1 in 5 infections recur²

1 in 7 open tibial fractures become infected^{3,4}

20 - 30% diabetic foot infections involve bone⁵

6 times higher cost when fracture-related infection is present⁶

Moreover, bone infection is associated with an increased risk of amputation. There is a call for effectively promoting and protecting bone healing when bone infection is present.

¹ Vasudevan, R., 'Biofilms : Microbial Cities of Scientific Significance', Journal of Microbiology and Experimentation, 1.3 (2014), 84–98.

² Conterno, L.O., Turchi, M.D. Antibiotics for treating chronic osteomyelitis in adults. Cochrane Database Syst Rev. 2013;2013(9).

³ Hoekstra, H., Smeets, B., Metsemakers, W-J., Spitz, A-C., Nijs, S. Economics of open tibial fractures: the pivotal role of length-of-stay and infection. Health Econ Rev. 2017 Sep;7(1):32.

⁴ Papakostidis, C., Kanakaris, N.K., Pretel, J., Faour, O., Morell, D.J., and Giannoudis, P.V., 'Prevalence of Complications of Open Tibial Shaft Fractures Stratified as per the Gustilo-Anderson Classification', Injury, 42.12 (2011), 1408–15.

⁵ Prompers, L., Huijberts, M., Apelqvist, J., Jude, E., Piaggese, A., Bakker, K., et al., 'High Prevalence of Ischaemia, Infection and Serious Comorbidity in Patients with Diabetic Foot Disease in Europe. Baseline Results from the Eurodiale Study', Diabetologia, 50.1 (2007), 18– 25.

⁶ O'Connor, O., Thahir, A., Krkovic, M., 'How Much Does an Infected Fracture Cost?', The Archives of Bone and Joint Surgery, 135.10(2) (2022), 135–40.




MANAGING BONE INFECTION WITH CERAMENT®

Promote bone remodeling

The traditional way of managing bone infection is to first clear the infection site from any dead and/or poorly perfused bone by surgical debridement, thus leaving a void in the bone.

CERAMENT® technology has the ideal characteristics for a bone substitute to manage bone voids and gaps:

- It resorbs at the same rate as bone naturally forms because of its unique composition (see diagram to the right)
- It is injectable into bone voids so that there is no dead space for bacteria to grow in

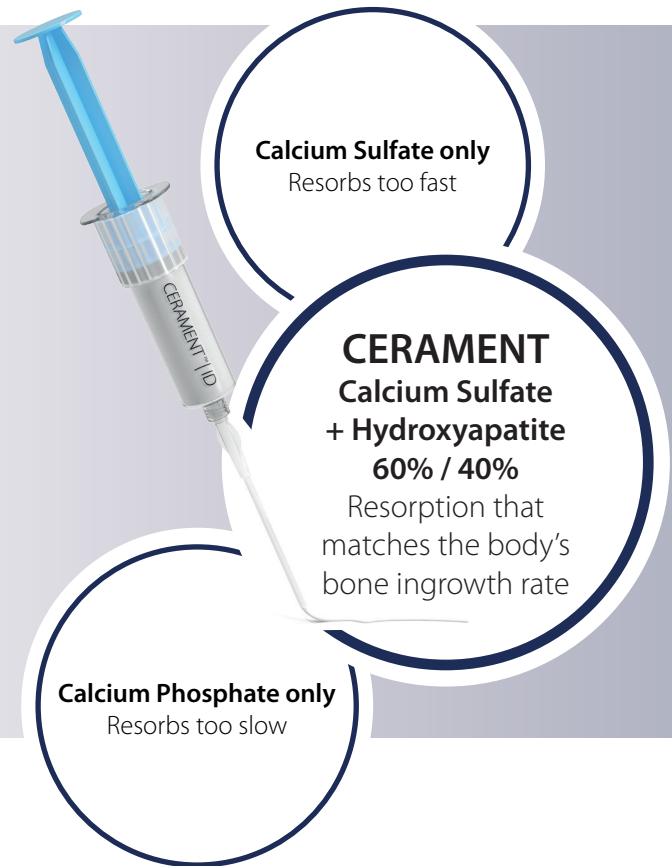


Protect bone healing

CERAMENT G and CERAMENT V are the antibiotic-eluting variants of the CERAMENT technology*. Our patented mixing system is specially designed to enable the safe and reliable local elution of antibiotics, gentamicin and vancomycin respectively.

Using CERAMENT G and CERAMENT V protects bone voids against colonization with a broad range of gram positive and gram negative (CERAMENT G only) bacteria.

*CERAMENT Bone Void Filler (BVF) is the third product of the CERAMENT technology which does not contain antibiotics

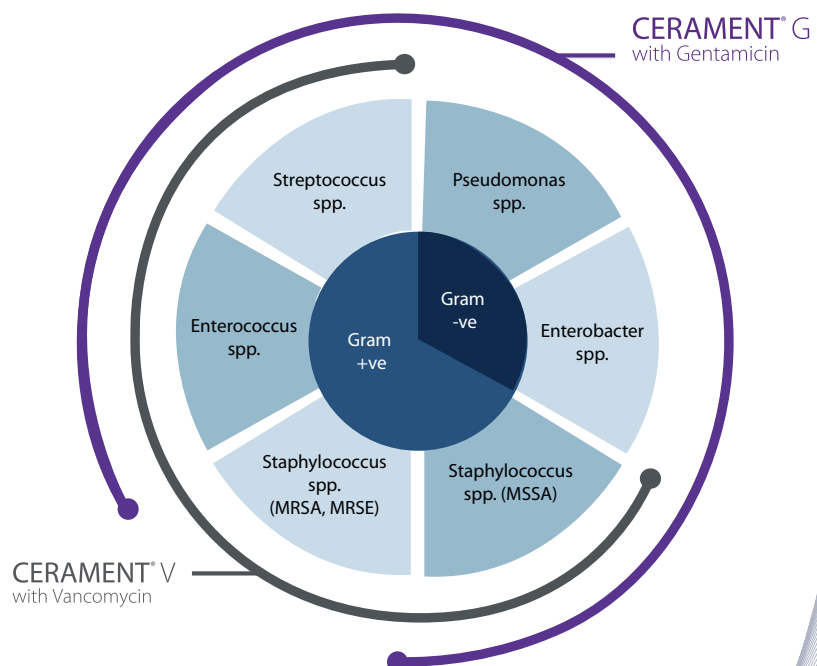


Gentamicin is a broad-spectrum antibiotic that is active against both gram positive and gram negative bacteria. Vancomycin is active against gram positive bacteria only. Spp. means 'species pluralis' or multiple species.

ALL IN ONE: CLOSED MIXING SYSTEM IN ONE COMPLETE KIT

How this system improves efficiency in operating theatres:

- All-in-one kit for bone void filling
- All surface sterile
- Standardized mixing procedure
- Reproducible elution of local antibiotic
- Customizable filling using tip extenders
- Overall less exposure to antibiotics to staff

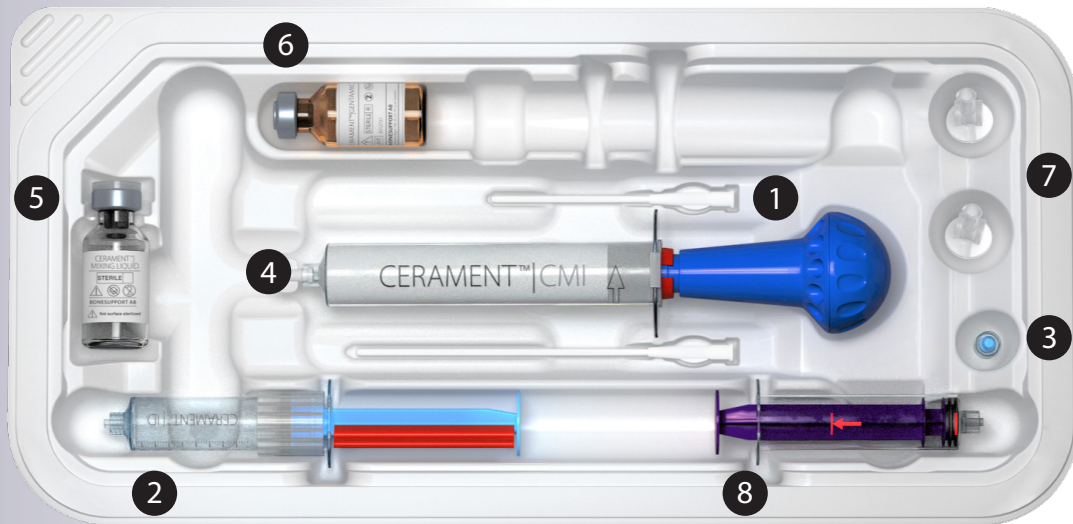

CERAMENT G kit

- | | |
|---------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|-------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| <ol style="list-style-type: none"> 1. 2 x 11G tip extenders with tapered ends in 50mm & 100mm lengths 2. Injection Device (ID) syringe 3. Valve 4. Combined Mixing and Injection (CMI) syringe pre-filled with hydroxyapatite (HA)/calcium sulfate (CaS) powder | <ol style="list-style-type: none"> 5. Mixing Liquid, sodium chloride 9 mg/mL liquid 6. Gentamicin sulfate, provides 17.5 mg gentamicin/mL paste 7. 2 x Dispensing pins 8. Syringe for preparing the gentamicin solution |
|---------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|-------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|

CERAMENT G



Size availability
5ml, 10ml


CERAMENT V



Size availability
10ml

No additional mixing equipment needed.

ROBUST EVIDENCE ON CERAMENT® TECHNOLOGY

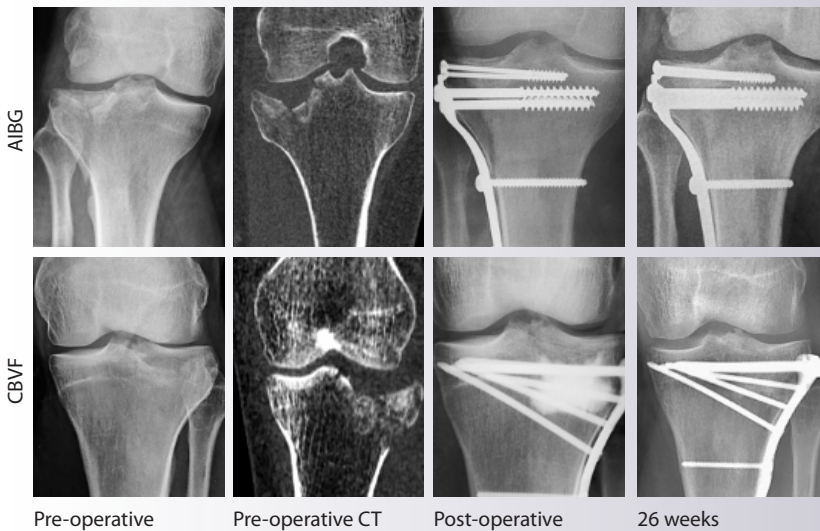

Evidence-based bone healing



- Largest amount of pre-clinical and clinical data of any bone substitute on the market (240+ publications and counting)
- Level 1 randomized controlled trial (RCT) – the CERTIFY study involves 135 patients and shows that CERAMENT Bone Void Filler (BVF) is as good as autograft (AIBG) in bone remodeling¹
- One and only orthobiologic product with robust long-term evidence: 94% infection-free after 6 year follow-up²



Safe and reliable antibiotic elution

- High burst of local antibiotic with sustained release above the minimum inhibitory concentration (MIC) for at least 28 days (see elution curves)^{3,4}
- No evidence of renal impairment caused by local gentamicin elution up to 525 mg⁵

CERAMENT G and CERAMENT V elution curves from *in vitro* data

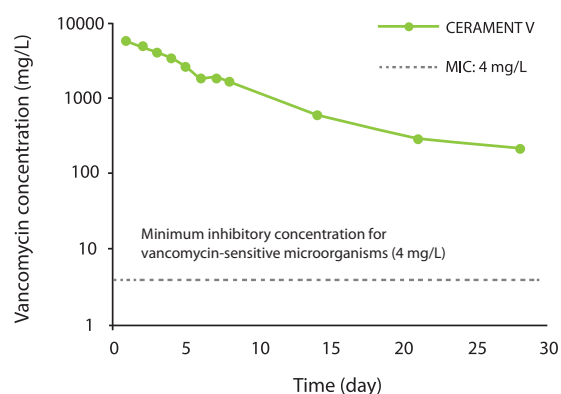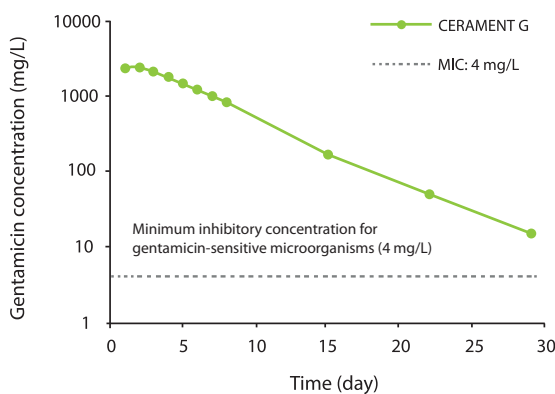

¹ Hofmann, A., Gorbulev, S., Guehring, T., Schulz, A.P., Schupfner, R., Raschke, M., et al., 'Autologous Iliac Bone Graft Compared with Biphasic Hydroxyapatite and Calcium Sulfate Cement for the Treatment of Bone Defects in Tibial Plateau Fractures: A Prospective, Randomized, Open-Label, Multicenter Study', The Journal of Bone and Joint Surgery. American Volume, 102.3 (2020), 179–93

² McNally, M.A., Ferguson, J.Y., Scarborough, M., Ramsden, A., Stubbs, D.A., Atkins, B.L., 'Mid- to Long-Term Results of Single-Stage Surgery for Patients with Chronic Osteomyelitis Using a Bioabsorbable Gentamicin-Loaded Ceramic Carrier', The Bone & Joint Journal, 104-B.9 (2022), 1095–1100

³ Stravinskas, M., Horstmann, P., Ferguson, J.Y., Hettwer, W., Nilsson, M., Tarasevicius, S., et al., 'Pharmacokinetics of Gentamicin Eluted from a Regenerating Bone Graft Substitute', Bone and Joint Research, 5.9 (2016), 427–35

⁴ CERAMENT G. Data on File (MC-002489). CERAMENT V. Data on File (S031/2013). BONESUPPORT AB Sweden.

⁵ Muir, R., Birnie, C., Hyder-Wilson, R., Ferguson, J., McNally, M.A., 'Does Local Implantation of Gentamicin Impair Renal Function in Patients Undergoing Surgery for Chronic Bone Infection?', International Journal of Research in Orthopaedics, 7.3 (2021), 438

WHAT IS VALUE?

Value is about achieving the best outcomes per money spent.

CERAMENT G and CERAMENT V provide great value by enabling strategies to reduce surgical interventions and to increase efficiency in clinical practice. Taking advantage of their unique set of features paves the way for more cost-effective and patient-centered care.


Unique features of CERAMENT G and CERAMENT V

Feature	Importance
Proven bone remodeling	Patients are at greater risk of fractures and (re)infection when there is insufficient bone growth in a void
Proven consistent antibiotic elution above minimum inhibitory concentration (MIC) for a clinically relevant time period	Prolonged low-level exposure to antibiotics may cause antibiotic-resistant organisms to evolve
All-in-one kit for closed mixing and injection (antibiotic included)	One kit used straight off the shelf Overall less exposure to antibiotics No off-label mixing
Reproducible mix & setting times	Standardized procedures minimize errors and ensures consistent delivery
Injectable to completely fill bone voids	Dead space left in a bone void are conducive environments for residual bacteria to grow, thus increasing the risk of infection recurrence
Evidence of sustained clinical effectiveness beyond six years	Bone infection can recur beyond two years of follow-up

IMPROVED CLINICAL OUTCOMES


Against PMMA beads + autograft

Injectable CERAMENT G (CG) paste vs. Gentamicin-impregnated PMMA* beads + autograft with comparable patient cohorts^{1, 2}

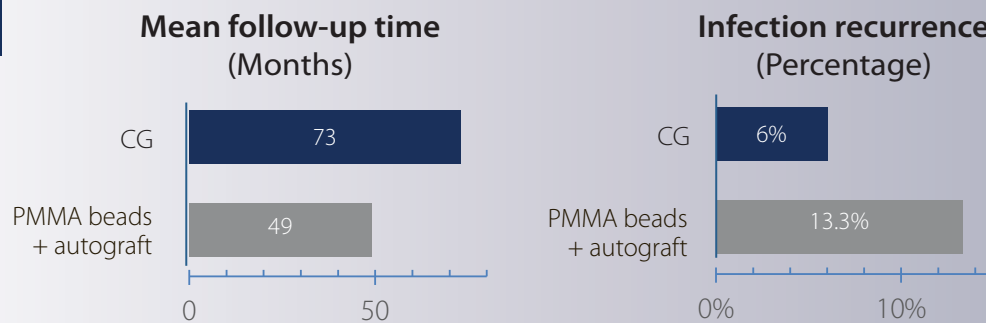

Against calcium sulfate (CS) pellets containing tobramycin

Injectable CERAMENT G (CG) paste vs. calcium sulfate (CS) pellets containing tobramycin³

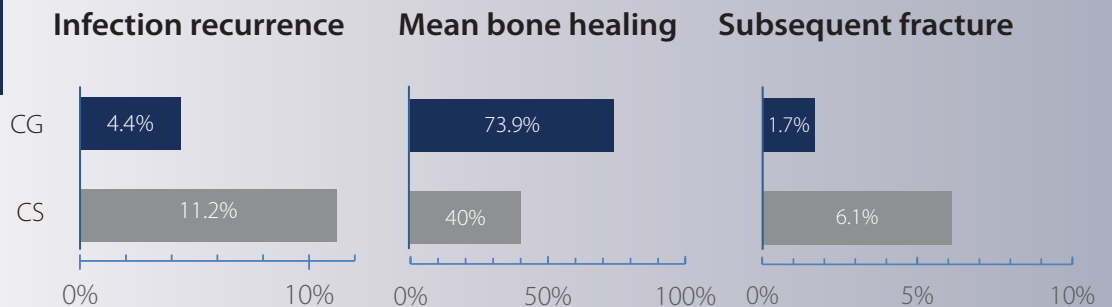

*Polymethyl methacrylate (PMMA)

Evidence of cost savings

A UK health economic study compares the healthcare utilization of a single-stage protocol using CERAMENT G at the Oxford bone infection unit versus the top 10 busiest osteomyelitis centers in England⁴.

17 less bed days
£8,500 / €10,043/CA\$14,883
direct cost savings per patient
from reduced length of stay*

*Based on an average cost per bed day of £500; the average annual exchange rate between pound sterling, euro and Canadian dollar for 2024 was used (£1 = €1.1633); data sourced from the Office for National Statistics (ons.gov.uk).

¹ McNally, M.A., Ferguson, J.Y., Scarborough, M., Ramsden, A., Stubbs, D.A., Atkins, B.L., 'Mid- to Long-Term Results of Single-Stage Surgery for Patients with Chronic Osteomyelitis Using a Bioabsorbable Gentamicin-Loaded Ceramic Carrier; The Bone & Joint Journal, 104-B.9 (2022), 1095–1100
² McNally, M. A., Small, J.O., Tofghi, H.G., Mollan R.A.B., 'Two-Stage Management of Chronic Osteomyelitis of the Long Bones; British Editorial Society of Bone and Joint Surgery, 75.3 (1993), 375–80
³ Ferguson, J., Bourget-Murray, J., Stubbs, D., McNally, M., Hotchen, A.J., 'A Comparison of Clinical and Radiological Outcomes between Two Different Biodegradable Local Antibiotic Carriers Used in the Single-stage Surgical Management of Long Bone Osteomyelitis; Bone & Joint Research, 12.7 (2023), 412–22
⁴ Ferguson, J., Alexander, M., Bruce, S., O'Connell, M., Beecroft, S., McNally, M.A., 'A Retrospective Cohort Study Comparing Clinical Outcomes and Healthcare Resource Utilisation in Patients Undergoing Surgery for Osteomyelitis in England: A Case for Reorganising Orthopaedic Infection Services; Journal of Bone and Joint Infection, 6.May (2021), 151–63

ENABLING SINGLE-STAGE PROTOCOL

Add convenience to your practice and remove the need for additional surgeries.

TRADITIONAL TREATMENT

Multi-stage surgery plus systemic antibiotics to manage infections

One complete kit with CERAMENT G or CERAMENT V

Single-stage surgery to manage infection

0-2

WEEKS

3-4

WEEKS

5-6

WEEKS

7+

WEEKS

MULTI STAGES

1st surgical procedure
Placement of non- absorbable antibiotic carrier.

HEALING STARTS

Subsequent procedures
Harvest bone graft (for autograft), removal of non-resorbable antibiotic carrier plus bone void filling with bone graft.

HEALING STARTS

Single-stage surgery plus systemic antibiotics to manage infections with CERAMENT G or CERAMENT V

Single stage protocols featuring CERAMENT G and CERAMENT V

Name	Clinic	Indication	Description
The Oxford Protocol^{1, 5}	Nuffield Orthopaedic Centre, Oxford University Hospitals, NHS, United Kingdom	Cavitary and segmental bone infection	Biphasic bioceramics such as CERAMENT G and CERAMENT V are assigned for cortico-medullary and combined defects. Evidence of cost savings with protocol
The Fix and Flap^{2, 3, 6}	Wythenshawe Hospital, Manchester University NHS Foundation Trust, United Kingdom	Open traumatic fractures	An ortho-plastic approach where definitive skeletal stabilization and soft tissue coverage are conducted in the same sitting. CERAMENT G is used as adjuvant local antibiotic bio-composite
The Silo Technique⁴	Wythenshawe Hospital, Manchester University NHS Foundation Trust, United Kingdom	Diabetic foot osteomyelitis	Debridement of the dead bone and local delivery of antibiotics in drilled tunnels using CERAMENT G

¹ Mifsud, M., McNally, M.A., 'Local Delivery of Antimicrobials in the Treatment of Bone Infections', Orthopaedics and Trauma, 33.3 (2019), 160-65

² Jahangir, N., Niazi, N., Aljawadi, A., Reid, A., Wong, J., Drampalos, E., et al., 'The Use of Adjuvant Local Antibiotic Hydroxyapatite Bio- Composite in the Management of Open Gustilo Anderson Type IIIB Fractures. A Prospective Review', Journal of Orthopaedics, 16.3 (2019), 278-82

³ Aljawadi, A., Islam, A., Jahangir, N., Niazi, N., Ferguson, Z., Sephton, B., et al., 'Adjuvant Local Antibiotic Hydroxyapatite Bio- Composite in the Management of Open Gustilo Anderson IIIB Fractures. Prospective Review of 80 Patients from the Manchester Ortho-Plastic Unit', Journal of Orthopaedics, 18.December 2019 (2020), 261-66


⁴ Drampalos, E., Mohammad, H.R., Kosmidis, C., Balal, M., Wong, J., Pillai, A., 'Single Stage Treatment of Diabetic Calcaneal Osteomyelitis with an Absorbable Gentamicin-Loaded Calcium Sulphate/Hydroxyapatite Biocomposite: The Silo Technique', Foot, 34 (2018), 40-44

⁵ Ferguson, J., Alexander, M., Bruce, S., O'Connell, M., Beecroft, S., McNally, M.A., 'A Retrospective Cohort Study Comparing Clinical Outcomes and Healthcare Resource Utilisation in Patients Undergoing Surgery for Osteomyelitis in England: A Case for Reorganising Orthopaedic Infection Services', Journal of Bone and Joint Infection, 6.May (2021), 151-63

⁶ Henry, J.A., Ali, A., Elkhidir, I.H., Reid, A., Wong, J., Pillai, A., 'Long-Term Follow-Up of Open Gustilo-Anderson IIIB Fractures Treated With an Adjuvant Local Antibiotic Hydroxyapatite Bio-Composite', Cureus, 15.Mic (2023)

ENABLING SINGLE-STAGE PROTOCOL


Example of filling a bone void with CERAMENT G in a single-stage protocol

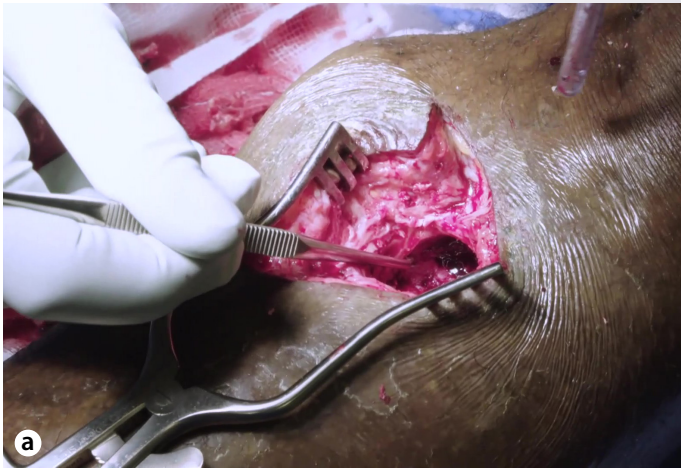

a Dead bone is removed leaving dead space.



b CERAMENT G is injected into the bone void.

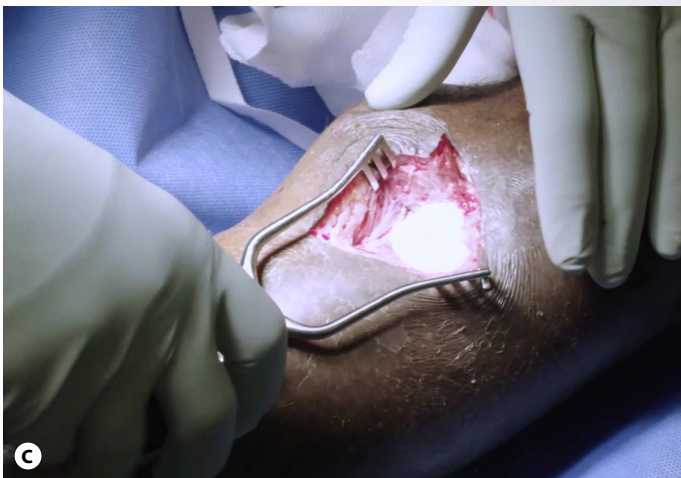


c Bone void is completely filled.



d CERAMENT G has set and wound is closed.

Note: Images are from a chronic osteomyelitis case and they are reproduced with kind permission from the Nuffield Orthopaedic Centre, Oxford UK.

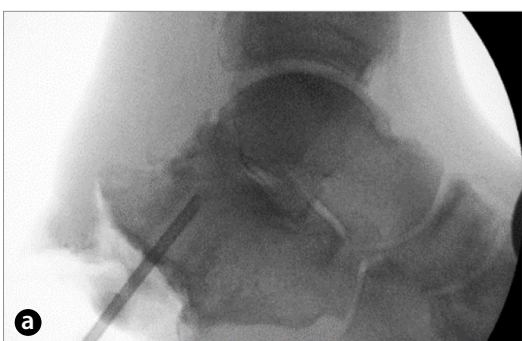

The Silo Technique.
Silo type tunnels (four to five) are drilled into the os calcis and filled with 5ml of CERAMENT G¹

¹ Drampalos, E., Hasan R.M., Kosmidis, C., Balal, M., Wong, J., Pillai, A., 'Single Stage Treatment of Diabetic Calcaneal Osteomyelitis with an Absorbable Gentamicin-Loaded Calcium Sulphate/Hydroxyapatite Biocomposite: The Silo Technique', Foot, 34 (2018), 40-44

OUTCOMES OF SINGLE-STAGE PROTOCOLS

Clinical outcomes of single-stage protocols that includes CERAMENT G and CERAMENT V

The Oxford Protocol¹
Chronic osteomyelitis
 100 patients

Mean follow-up:
 6.05 years

94%
 remained infection-free

97%
 did not develop
 a fracture

The Fix and Flap²
Fracture-related infections
 81 patients

Mean follow-up time:
 55.8 months

96.3%
 deep infection free

96.3%
 limb salvage rate

96%
 bony union rate

Retrospective study comparing local antibiotic using CERAMENT G or CERAMENT V therapy vs current standard of care without local antibiotic therapy³

Diabetic foot osteomyelitis
 56 patients treated with local antibiotic
 49 received standard of care without local antibiotic

Mean follow-up time:
 14 months

87,5%
 Survival rate

98%
 did not have a
 major amputation

CERAMENT G and CERAMENT V are valuable adjunctive therapies to keep a low infection recurrence rate in bone infection cases.

Potential cost savings from preventing (re)infection

Location of fracture	Savings = average cost of infected fracture ⁴ - average cost of uninfected fracture ⁴		
	Euros (€)	Pounds (£)	Canadian dollars (CA\$)
Tibia	35,758	30,265	52,991
Femur	20,322	17,200	30,115
Hip	16,278	13,777	24,122

Additionally, bone infection is associated with an increased risk of amputation. For patients with diabetic foot, amputation-related costs can range from £23,100/€25,550/CA\$41,300 to £29,700/€32,850/CA\$53,100^{5*}.

*Adjusted by PPP exchange rates for 2025 (US\$1=€0.66, US\$1=€0.73, US\$1=CA\$1.18); data sourced from the International Monetary Fund (imf.org)

¹ McNally, M.A., Ferguson, J.Y., Scarborough, M., Ramsden, A., Stubbs, D.A., Atkins, B.L., 'Mid- to Long-Term Results of Single-Stage Surgery for Patients with Chronic Osteomyelitis Using a Bioabsorbable Gentamicin-Loaded Ceramic Carrier', *The Bone & Joint Journal*, 104-B.9 (2022), 1095–1100

² Henry, J.A., Ali, A., Elkhidir, I.H., Reid, A., Wong, J., Pillai, A., 'Long-Term Follow-Up of Open Gustilo-Anderson IIIB Fractures Treated With an Adjuvant Local Antibiotic Hydroxyapatite Bio-Composite', *Cureus*, 15.Mic (2023)

³ Metaoy, S., Rusu, I., Pillai, A., 'Adjuvant Local Antibiotic Therapy in the Management of Diabetic Foot Osteomyelitis', *Clinical Diabetes and Endocrinology*, 2024, 1–6

⁴ O'Connor, O., Thahir, A., Krkovic, M., 'How Much Does an Infected Fracture Cost?', *The Archives of Bone and Joint Surgery*, 135.10(2) (2022), 135–40.

⁵ Petrakis, I., Kyriopoulos, I.J., Ginis, A., Athanasakis, K., 'Losing a Foot versus Losing a Dollar; a Systematic Review of Cost Studies in Diabetic Foot Complications', *Expert Review of Pharmacoeconomics and Outcomes Research*, 17.2 (2017), 165–80

A COST-EFFECTIVE STRATEGY FOR HOSPITALS

Get it right the first time - enabling single-stage protocols with CERAMENT G and CERAMENT V reduces the number of surgical procedures and has shown to improve clinical outcomes beyond six years

Having low infection rates can generate cost savings – evidence of 17 less bed days on average following the Oxford Protocol with CERAMENT G (£8,500/€10,043/CA\$14,883 saving)

CERAMENT G and CERAMENT V are all-in-one closed mixing systems – standardized mixing procedure minimizes error and ensures more consistency in bone void filling, overall less exposure to antibiotics to staff

BONESUPPORT is a Swedish orthobiologics company founded in 1999 that develops, manufactures, and markets CERAMENT® – an innovative portfolio of injectable and antibiotic-eluting bone void fillers.

Availability of CERAMENT® is dependent on its regulatory status in individual markets, contact your local representative.

For complete product information, including indications, contraindications, warnings, precautions and potential adverse events, see package insert.



CERAMENT® G with Gentamicin

A synthetic bone substitute with hydroxyapatite, calcium sulfate and the antibiotic gentamicin. Remodels into bone within 6-12 months and elutes gentamicin.



CERAMENT® V with Vancomycin

A synthetic bone substitute with hydroxyapatite, calcium sulfate and the antibiotic vancomycin. Remodels into bone within 6-12 months and elutes vancomycin.

Ordering details: E: order@bonesupport.com T: +46 (0)46 286 53 70